




eva
N E X U S



Remzi Avci

MD, Professor
Head of Retina Eye Hospital Bursa,
Turkey

"One of the features of EVA NEXUS that I find most useful in my clinical practice is its infusion reversal function. This means that extra pressure applied to the eye, for example during dyes, PFCL or any other liquid injection or indentation, will be released by pump reversal. This system is very efficient with the combination of the AVETA trocar cannula system, which maximises flow rates thanks to an uninterrupted flow path between the infusion line and the trocar cannula. The very fast VTi pump can increase and release pressure very fast due to small chambers driven by pistons. Whenever needed, the infusion part of the EVA NEXUS pump can switch to reverse mode, displacing liquid from the eye back to infusion line, to release excess pressure. This system helps you keep stable IOP during the full surgery regardless of the aspiration or injection procedures."

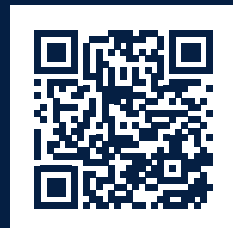

More info:
dorcglob.com/eva-nexus